Neu ab:

New as of:

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10.2018

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10.2018

CEREC camera

Care, cleaning, disinfection and sterilization

Instructions for CEREC Bluecam and CEREC Omnicam
(not valid for USA/Canada)

English

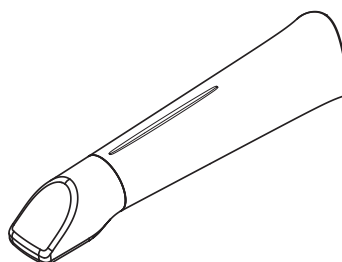
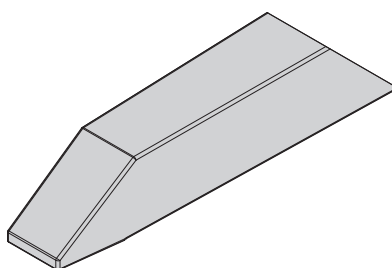
Table of contents

1	Cleaning agents and disinfectants	3
2	CEREC Omnicam	4
2.1	General information	4
2.2	Components of the Omnicam	5
2.3	Removing the mirror sleeve	5
2.4	Refitting the mirror sleeve	5
2.5	Pre-cleaning of the mirror sleeve	6
2.6	Wipe disinfection for the camera and mirror sleeve	6
2.7	High-level disinfection of the mirror sleeve with the HLD set 66 05 120	7
2.8	High-level disinfection of the mirror sleeve with the HLD set 63 46 907	11
2.9	Hot air sterilization of the mirror sleeve	13
2.10	Use of CEREC Omnicam disposable sleeves	14
3	CEREC Bluecam	15
3.1	CEREC Bluecam without sapphire glass	15
3.1.1	General information	15
3.1.2	Components of the CEREC Bluecam without sapphire glass	16
3.1.3	Removing the prismatic tube	16
3.1.4	Refitting the prismatic tube	17
3.1.5	Pre-cleaning of the prismatic tube	17
3.1.6	Wipe disinfection for the camera and prismatic tube	17
3.1.7	High-level disinfection of the prismatic tube without sapphire glass	18
3.1.8	Hot air sterilization of the prismatic tube	20
3.2	CEREC Bluecam with sapphire glass	21
3.2.1	General information	21
3.2.2	Components of the CEREC Bluecam with sapphire glass	21
3.2.3	Removing the prismatic tube	22
3.2.4	Refitting the prismatic tube	22
3.2.5	Pre-cleaning of the prismatic tube	22
3.2.6	Wipe disinfection for the camera and prismatic tube	23
3.2.7	High-level disinfection of the prismatic tube with sapphire glass ...	24

1 Cleaning agents and disinfectants

NOTE

Approved care, cleaning, and disinfecting agents

Use only care, cleaning, and disinfecting agents approved by Dentsply Sirona!

Alpro	• Minuten Spray classic • Minuten Wipes • Plasti Sept • Plasti Sept Wipes
Merz	• Pursept-A
Dürr	• FD 312

2 CEREC Omnicam

CAUTION

If the CEREC Omnicam accidentally falls down, check to make sure that the camera window and coated sapphire glass are not damaged. If the CEREC Omnicam has been damaged, it must no longer be used on patients.

The CEREC Omnicam must be recalibrated.

NOTE

Do not under any circumstances sterilize the CEREC Omnicam or the camera cable.

NOTE

The mirror sleeves cannot be sterilized in the autoclave.

NOTE

For markets where the RKI* guidelines are to be observed

The mirror sleeve is classified as a "semicritical medical device A" according to RKI guidelines and therefore does not have to be autoclavable.

*RKI=Robert Koch Institute, Berlin (Germany).

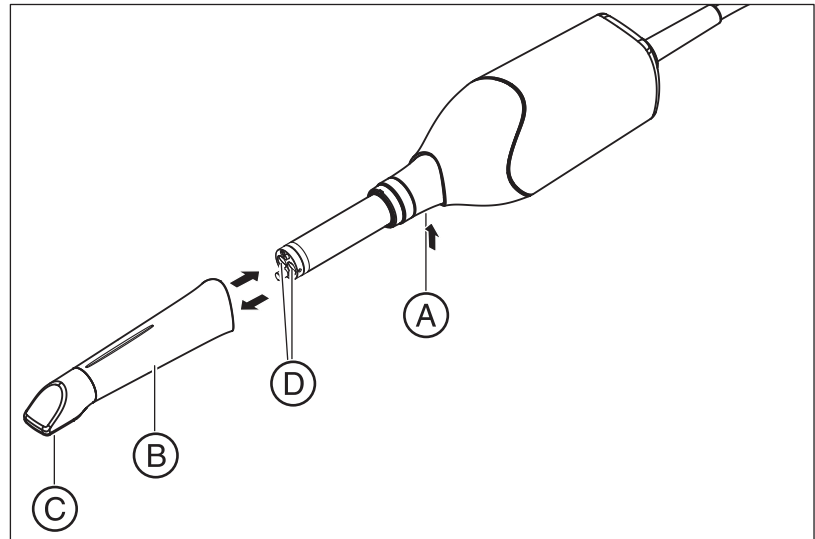
2.1 General information

The CEREC Omnicam is a very sensitive optical device and must therefore be handled with the **utmost care**. Protect the coated sapphire glass and the camera window against scratches and clean them with a lint-free cloth and ethanol (commercially available cleaning alcohol) if any haze is noted during the acquisition.

The following methods are available for cleaning the mirror sleeves:

- Wipe disinfection for the camera and mirror sleeve [→ 6]
- High-level disinfection of the mirror sleeve with the HLD set 66 05 120 [→ 7]
- High-level disinfection of the mirror sleeve with the HLD set 63 46 907 [→ 11]
- Hot air sterilization of the mirror sleeve [→ 13]

2.2 Components of the Omnicam

A	Press detent to release	C	Sapphire glass (coated)
B	Mirror sleeve	D	Camera window

2.3 Removing the mirror sleeve

If required to remove the mirror sleeve do the following:

1. Press the mirror sleeve against the camera body.
2. Press detent A.

NOTE

There is a risk of damaging the camera window or the coated sapphire glass if the mirror sleeve is not pushed straight toward the front.

> Push the mirror sleeve straight toward the front; **do not tilt it.**

3. Pull off the mirror sleeve.

2.4 Refitting the mirror sleeve

- > Carefully refit the mirror sleeve until it locks in place.

NOTE

There is a risk of damaging the camera window or the coated sapphire glass if the mirror sleeve is not pushed straight toward the front.

- > The mirror sleeve must not come into contact with the camera window.
- > Push the mirror sleeve straight toward the camera body; **do not tilt it.**

- > Carefully refit the mirror sleeve until it locks in place.

2.5 Pre-cleaning of the mirror sleeve

Clean the CEREC Omnicam immediately after use as follows:

1. Carefully wipe the mirror sleeve while it is on the camera so that no dirt whatsoever can remain stuck and can harden on the surface of the mirror sleeve. Use an absorbent cotton-based gauze for this with water and neutral soap or a lint-free cloth that has been soaked in one of the cleaning agents listed under "Cleaning agents and disinfectants [→ 3]".
2. Wipe down the mirror sleeve afterwards with the absorbent cotton gauze dipped in drinking water.
3. Next dry the mirror sleeve using a lint-free cloth.

2.6 Wipe disinfection for the camera and mirror sleeve

NOTE

Do not spray the CEREC Omnicam or immerse it in cleaning agents or disinfectants!

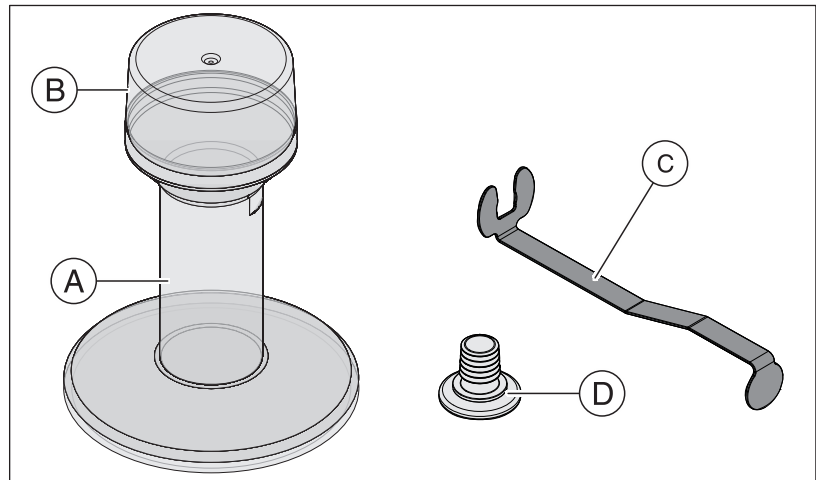
1. Carry out a preliminary cleaning process (see "Pre-cleaning of the mirror sleeve [→ 6]").
2. Use a new absorbent cotton gauze which is moistened with one of the agents listed in the "Cleaning agents and disinfectants [→ 3]" section. First wipe the camera enclosure and then the mirror sleeve. The camera enclosure and mirror sleeve are disinfected through this step.
3. Wipe down the mirror sleeve afterwards again with the absorbent cotton gauze dipped in drinking water.
4. Next dry the mirror sleeve using a lint-free cloth.

2.7 High-level disinfection of the mirror sleeve with the HLD set 66 05 120

NOTE

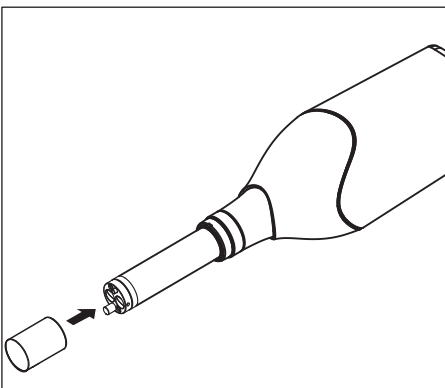
The hot air sterilization and high-level disinfection must not be combined.

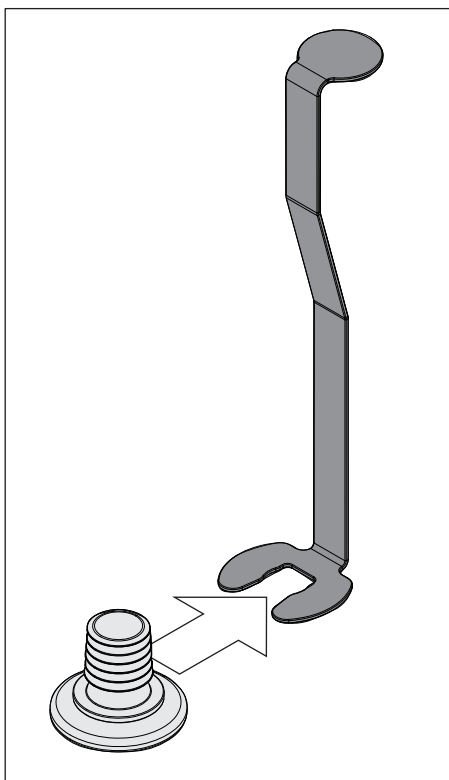
The complete process for high-level disinfection (HLD) is as follows – assuming the CIDEX[®] OPA disinfectant is available in your country via your dealer:

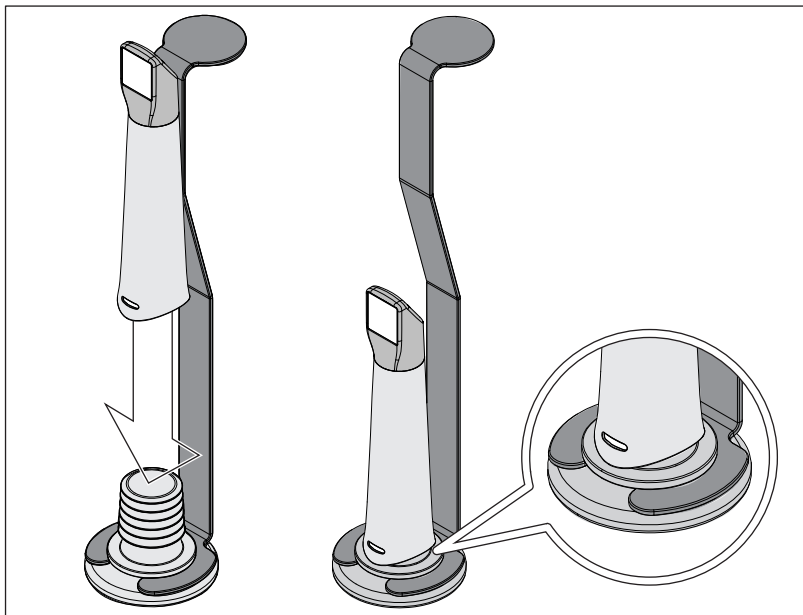
A	HLD tank	C	Metal bracket
B	Cover	D	Plugs

1. Carry out a preliminary cleaning process (see "Pre-cleaning of the mirror sleeve [→ 6]").
2. Always use personal protective equipment when implementing the high-level disinfection if not before.
3. Remove the mirror sleeve from the CEREC Omnicam (see "Removing the mirror sleeve [→ 5]").
4. Place the white protective cap on the camera head and place the camera in the camera storage unit.
5. Use the following disinfectant for the high-level disinfection: CIDEX[®] OPA.

6. Push the plug completely into the metal bracket.

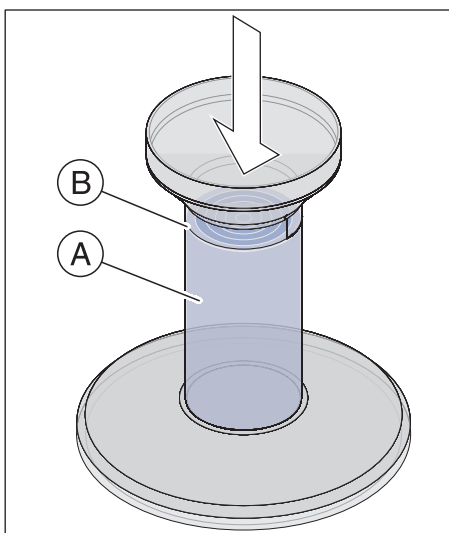
7. Press the mirror sleeve onto the plug while holding the metal bracket firmly in place in order to ensure that no particle contamination or fluids penetrate the inside of the mirror sleeve.

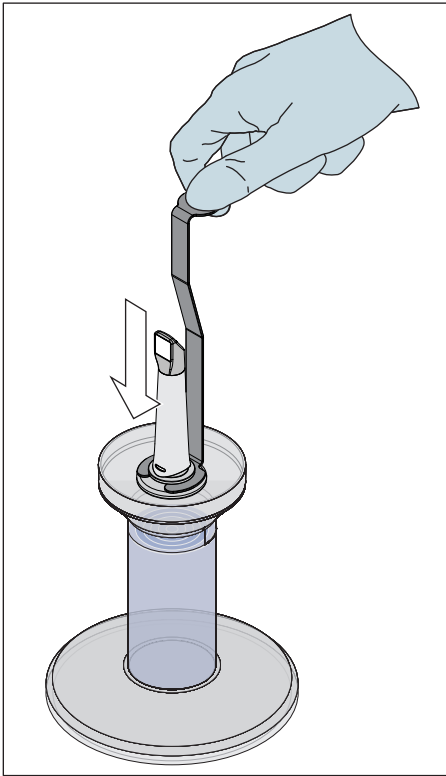
✎ The plug provides a waterproof seal for the mirror sleeve when attached correctly.

⚠ CAUTION

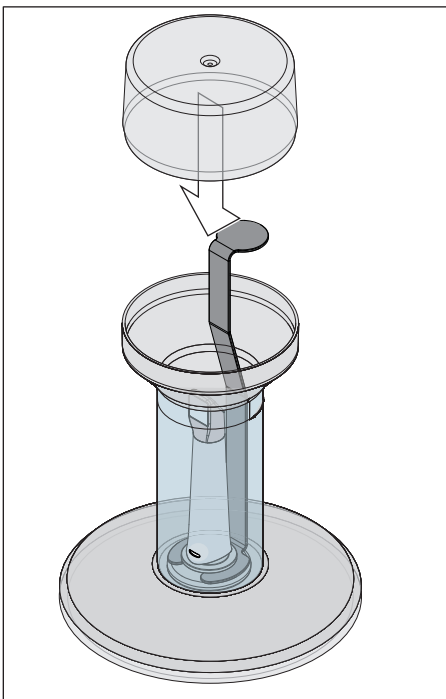
Observe the safety instructions from the disinfectant manufacturer.

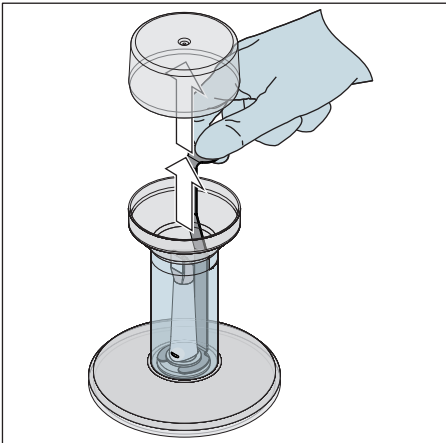
8. Carefully fill the HLD container up to the marking (B). You can use a funnel to pour in the fluid (A). Do not spill any disinfectant. Observe the safety instructions from the disinfectant manufacturer when cleaning if you do spill any disinfectant.

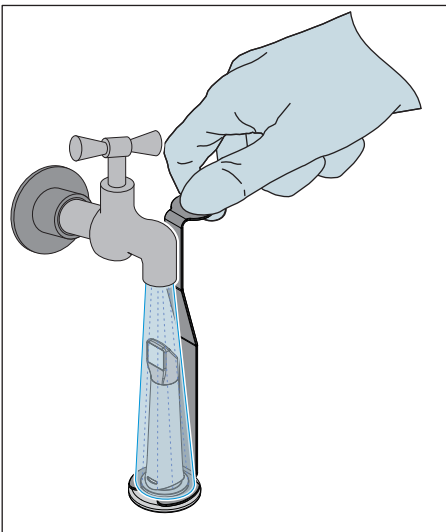
9. Insert the metal bracket with the metal sleeve.

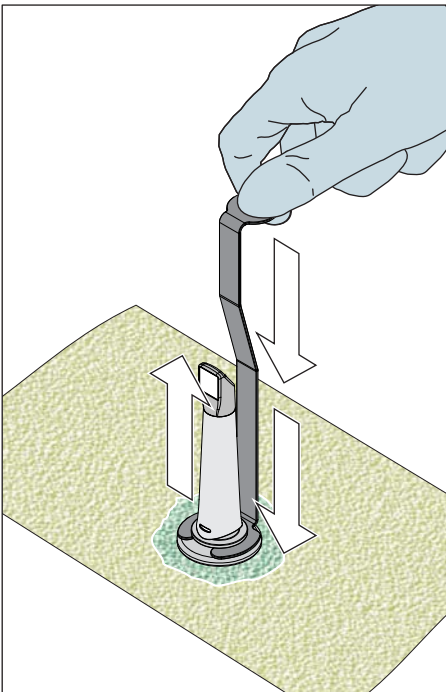
10. Place the cover on the container and leave the mirror sleeve in the disinfectant for at least 12 minutes (CIDEX® OPA). The mirror sleeve will not be damaged if it remains in the disinfectant for longer than the required time, but removing it shortly after the required time is recommended.

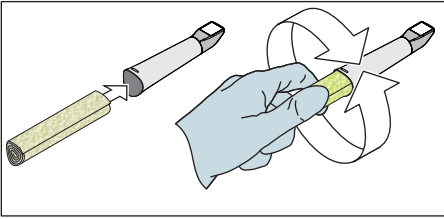
- 11.** Remove the cover and take the bracket with the mirror sleeve out of the container.

- 12.** Carefully rinse the mirror sleeve with tap water for at least 30 seconds while holding it with the bracket.

- 13.** Remove the bracket from the mirror sleeve carefully and slowly by pushing the bracket downwards. Ensure that the mirror sleeve is pointing upwards so that no fluids are able to penetrate into the interior of the mirror sleeve.
- 14.** Remove the plug from the bracket.
- 15.** Dry the plug.
- 16.** Dispose of the fluid and store the bracket in the empty container if the HLD set is not used for more than one week.

17. Dry the entire outer surface of the mirror sleeve as well as the inner section using a soft lint-free cloth. Ensure that no fluid is able to penetrate the inside of the mirror sleeve.
18. Store the mirror sleeve in such a way that it is protected from contamination until the next use.
19. Remove the white protective cap from the camera head before use.
20. Carefully re-attach the mirror sleeve and allow it to lock in place (see "Refitting the mirror sleeve [→ 5]").

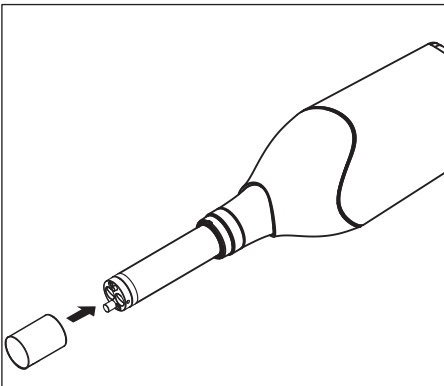
2.8 High-level disinfection of the mirror sleeve with the HLD set 63 46 907

NOTE

The hot air sterilization and high-level disinfection must not be combined.

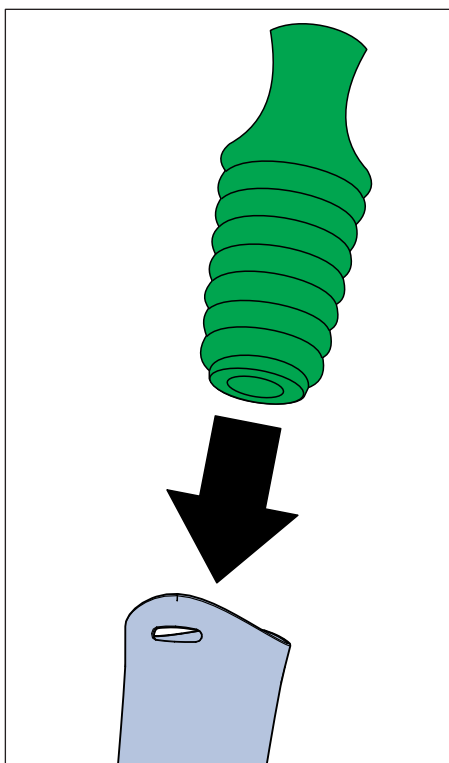
The complete process for high-level disinfection (HLD) is as follows – assuming the CIDEX[®] OPA disinfectant is available in your country via your dealer:

1. Carry out a preliminary cleaning process (see "Pre-cleaning of the mirror sleeve [→ 6]").
2. Always use personal protective equipment when implementing the high-level disinfection if not before.
3. Remove the mirror sleeve from the CEREC Omnicam (see "Removing the mirror sleeve [→ 5]").
4. Place the white protective cap on the camera head and place the camera in the camera storage unit.
5. Use the following disinfectant for the high-level disinfection: CIDEX[®] OPA.

⚠ CAUTION

Observe the safety instructions from the disinfectant manufacturer.

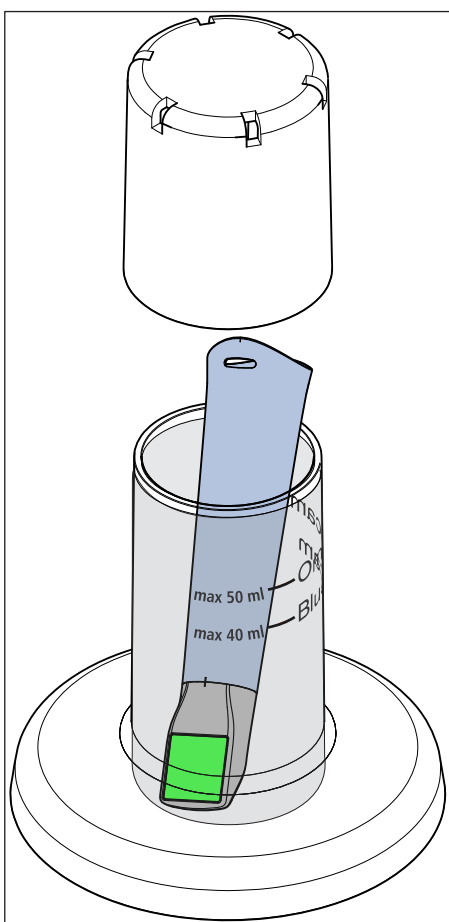
6. Press the green plug into the interior of the mirror sleeve in order to ensure that no particle contamination or fluids penetrate the inside of the mirror sleeve.

NOTE

Nature and source of the danger

The plug does not seal the mirror sleeve watertight.

7. Carefully fill the HLD container up to the 50 ml mark ("Omnicam" mark). A funnel can be used to help pour in the liquid. Do not spill the disinfectant. Should you spill any disinfectant, follow the safety instructions from the disinfectant manufacturer.

NOTE

Nature and source of the danger

Only place the mirror sleeve into the container and not the entire camera.

Ensure that the mirror sleeve is placed upright in the disinfectant (with the coated sapphire glass down) so that the solution cannot penetrate the inside of the mirror sleeve.

8. Insert the mirror sleeve in the correct position (coated sapphire glass down).
9. Place the cover on the container and leave the mirror sleeve in the disinfectant for at least 12 minutes (CIDEX[®] OPA). The mirror sleeve will not be damaged if it remains in the disinfectant for longer than the required time, but removing it shortly after the required time is recommended.
10. Remove the mirror sleeve from the container.
11. Carefully rinse the mirror sleeve with tap water for at least 30 seconds. Ensure that no water is able to penetrate into the interior of the mirror sleeve.
12. Wipe down the mirror sleeve. Use cotton wool gauze to do so, moistened with tap water.
13. Dry the mirror sleeve using a soft, lint-free cloth.
14. Remove the green plug from the mirror sleeve.

15. Store the mirror sleeve in such a way that it is protected from contamination until the next use.
16. Remove the white protective cap from the camera head before use.
17. Carefully re-attach the mirror sleeve and allow it to lock in place (see "Refitting the mirror sleeve [→ 5]").

2.9 Hot air sterilization of the mirror sleeve

NOTE

Hot air sterilization and high-level disinfection must not be combined.

The process for hot air sterilization is as follows:

1. Carry out a preliminary cleaning process (see "Pre-cleaning of the mirror sleeve [→ 6]").
2. Remove the mirror sleeve from the CEREC Omnicam (see "Removing the mirror sleeve [→ 5]" section).
3. The mirror sleeve can be sterilized using hot air (180 °C, 30 min). Please ask your dealer for deals on hot air sterilizers. Place the mirror sleeves in the hot air sterilizer and follow the manufacturer's instructions.
4. Store the mirror sleeve in such a way that it is protected from contamination until the next use.
5. Carefully re-attach the mirror sleeve and allow it to lock in place (see "Refitting the mirror sleeve [→ 5]" section).

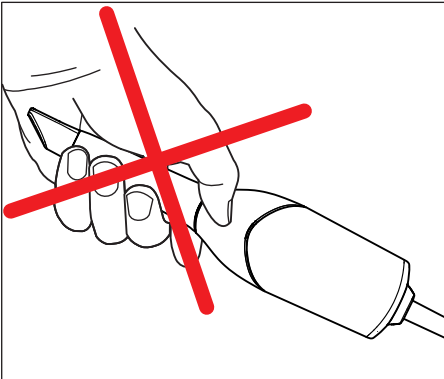
2.10 Use of CEREC Omnicam disposable sleeves

The CEREC Omnicam can be operated with disposable plastic sleeves in order to ensure maximum infection control. The sleeves are available using order number 66 32 264.

NOTE

The disposable sleeves cannot be used in the following circumstances:

- with the CEREC AF and CEREC AI system,
- when using color analysis,
- when using CEREC Ortho SW software.

1. Remove the disposable sleeve from the CEREC Omnicam reflective sleeve after each patient and dispose of this following standard procedures.
2. Disinfect the camera with a wipe (see "Wipe disinfection for the camera and mirror sleeve [→ 6]").
3. Remove one disposable sleeve from its packaging. Slide the sleeve onto the CEREC Omnicam reflective sleeve until it stops so that it is firmly held on the camera. Hold the camera in other places, not just on the sleeve in order to prevent it from falling.
4. Place the camera in the tray so that it can warm up for at least 15 minutes before intraoral use. Ensure that the tray is disinfected (see "Non-critical surfaces with no monitor").
5. When using the intraoral camera with disposable sleeves, the camera should be positioned closer to the teeth than it would be without a sleeve. If the 3D preview indicates artifacts, check whether the window of the disposable sleeve is clean and wipe it down with alcohol.

3 CEREC Bluecam

CAUTION

If the Bluecam has accidentally fallen down, check to make sure that the front lens and prism are not damaged. If it has been damaged, the CEREC Bluecam must not be used on patients any more.
The CEREC Bluecam must be recalibrated in any case.

NOTE

Do not under any circumstances sterilize the CEREC Bluecam or the video cable!

NOTE

The prism tubing is not suitable for autoclaving!

NOTE

For markets where the RKI* guidelines are to be observed

The prismatic tube is classified as a "semicritical medical device A" according to RKI guidelines and therefore does not have to be autoclavable.

*RKI=Robert Koch Institute, Berlin (Germany).

3.1 CEREC Bluecam without sapphire glass

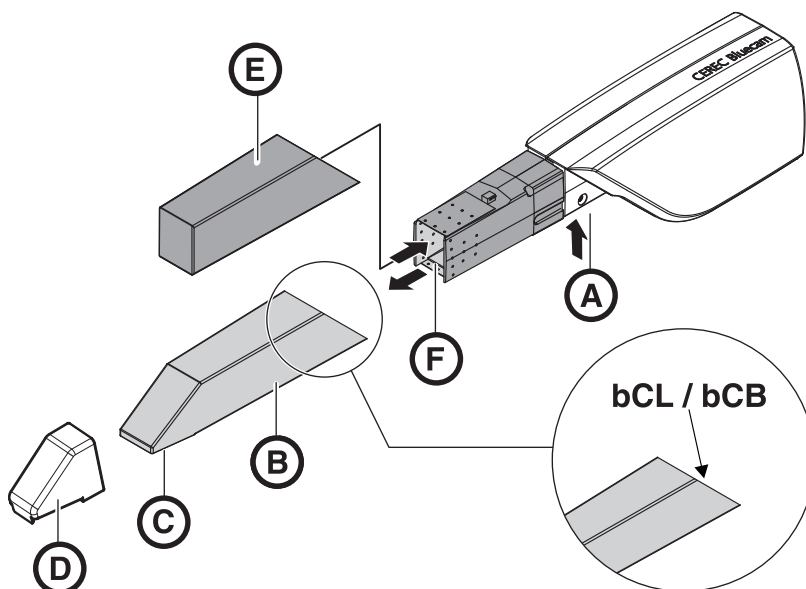
3.1.1 General information

The CEREC Bluecam is a very sensitive optical device and must therefore be handled with the **utmost care**. Protect the front lens or prism against scratches and clean them with a lint-free cloth and ethanol (commercially available cleaning alcohol) if any haze is noted during the acquisition.

The following methods are available for cleaning the prismatic tube:

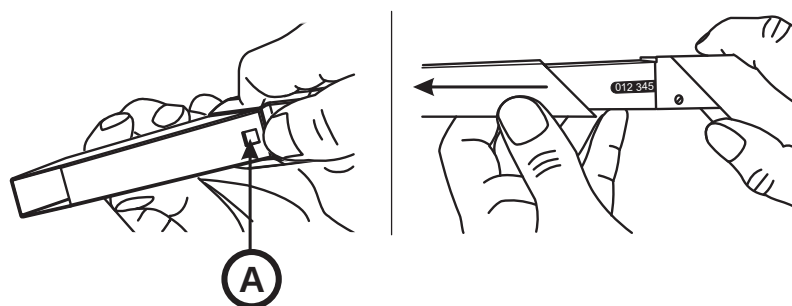
- Wipe disinfection for the camera and prismatic tube [→ 17]
- High-level disinfection of the prismatic tube without sapphire glass [→ 18]
- Hot air sterilization of the prismatic tube [→ 20]

3.1.2 Components of the CEREC Bluecam without sapphire glass

A	Press detent to release	D	Camera support 6 pcs, Order No. 5945360
B	Prismatic tube	E	Protective cap
C	Prism without sapphire glass (bCL / bCB)	F	Front lens

3.1.3 Removing the prismatic tube

A	Detent
---	--------

1. Press the prismatic tube against the camera body.
2. Press detent A.

NOTE

Risk of damaging the front lens or prism.

> Push the prismatic tube straight toward the front, do not tilt it.

3. Pull off the prismatic tube.

3.1.4 Refitting the prismatic tube

NOTE

Do not use CEREC 2 / CEREC 3 prismatic tubes.

NOTE

Risk of damaging the front lens or prism.

- > The prismatic tube must not touch the front lens.
- > Push the prismatic tube straight toward the camera body; do not tilt it.

- > Carefully refit the prismatic tube until it locks in place.

3.1.5 Pre-cleaning of the prismatic tube

Cleaning process immediately after using the CEREC Bluecam: Wipe the camera with the appropriate prismatic tube so that the superficial impurities cannot harden and stick to the surface of the prismatic tube. Remove any visible contaminations using cotton wool gauze, moistened with tap water and a neutral soap, or lint-free cloth soaked with one of the disinfectant agents specified in the "Cleaning agents and disinfectants [→ 3]" section. Wipe down the prismatic tube afterwards again with the absorbent cotton gauze dipped in drinking water. Next dry the prismatic tube using a lint-free cloth.

3.1.6 Wipe disinfection for the camera and prismatic tube

NOTE

Do not spray the CEREC Bluecam with or immerse it in cleaning agents or disinfectants!

1. Remove the camera support (if present).
2. Carry out a preliminary cleaning process (see "Pre-cleaning of the prismatic tube [→ 17]").
3. Use a new absorbent cotton gauze which is moistened with one of the agents listed in the "Cleaning agents and disinfectants [→ 3]" section. First wipe the camera enclosure and then the prismatic tube. The camera enclosure and prismatic tube are disinfected through this step.
4. Wipe down the prismatic tube afterwards again with the absorbent cotton gauze dipped in drinking water.
5. Next dry the prismatic tube using a lint-free cloth.

3.1.7 High-level disinfection of the prismatic tube without sapphire glass

NOTE

Hot air sterilization and high-level disinfection must not be combined.

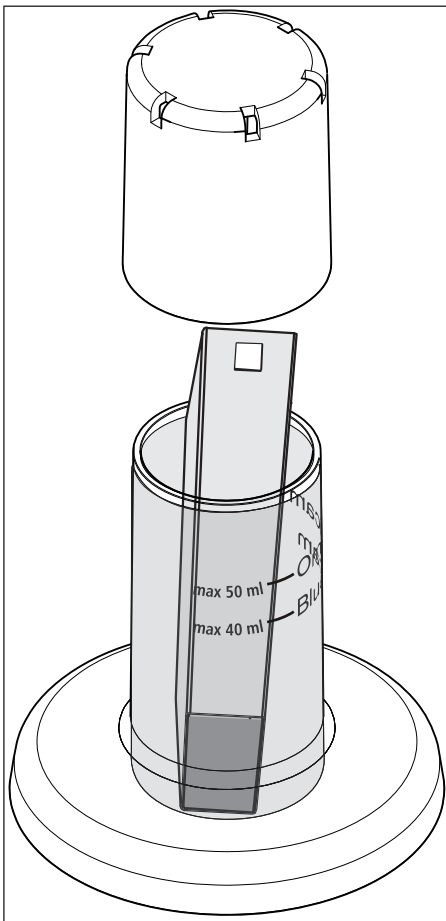
The complete process for high-level disinfection (HLD) is as follows – assuming the CIDEX[®] OPA disinfectant is available in your country via your dealer:

A HLD set to support the HLD processes can be ordered from Dentsply Sirona with REF 63 46 907.

1. Remove the camera support (if present).
2. Carry out a preliminary cleaning process (see "Pre-cleaning of the prismatic tube [→ 17]").
3. Always use personal protective equipment when implementing the high-level disinfection if not before.
4. Remove the prismatic tube from the CEREC Bluecam (see "Removing the prismatic tube [→ 16]").
5. Place the protective cap on the camera head and place the camera in the camera storage unit.
6. Use the following disinfectants for the high-level disinfection: CIDEX[®] OPA.

CAUTION

Observe the safety instructions from the disinfectant manufacturer.

7. Carefully fill the HLD container up to the 40 ml mark ("Bluecam" mark). A funnel can be used to help pour in the liquid. Do not spill the disinfectant. Should you spill any disinfectant, follow the safety instructions from the disinfectant manufacturer.

NOTE

Only place the prismatic tube into the container and not the entire camera.

Ensure that the prismatic tube is placed upright in the disinfectant (with the prism facing down) so that the solution cannot penetrate the inside of the prismatic tube.

8. Insert the prismatic tube in the correct position (with the prism facing down).
9. Place the cover on the container and leave the prismatic tube in the disinfectant for at least 12 minutes (CIDEX® OPA). The prismatic tube will not be damaged if it remains in the disinfectant for longer than the required time, but removing the prismatic tube shortly after the required time is recommended.
10. Remove the prismatic tube from the container.
11. Carefully rinse the prismatic tube with tap water for at least 30 seconds. Ensure that no water is able to penetrate into the interior of the prismatic tube.
12. Wipe down the prismatic tube. Use cotton wool gauze to do so, moistened with tap water.
13. Dry the prismatic tube using a soft, lint-free cloth.
14. Store the prismatic tube in such a way that it is protected from contamination until the next use.
15. Remove the protective cap from the camera head before use.
16. Carefully re-attach the prismatic tube and allow it to lock in place. Refer to the section on "Refitting the prismatic tube [→ 17]".

3.1.8 Hot air sterilization of the prismatic tube

NOTE

Hot air sterilization and high-level disinfection must not be combined.

The process for hot air sterilization is as follows:

1. Remove the camera support (if present).
2. Carry out a preliminary cleaning process (see "Pre-cleaning of the prismatic tube [→ 17]").
3. Remove the prismatic tube from the CEREC Bluecam (see the section on "Removing the prismatic tube [→ 16]").
4. The prismatic tube can be sterilized using hot air (180 °C, 30min). Please ask your dealer for deals on hot air sterilizers. Place the prismatic tube in the hot air sterilizer and follow the manufacturer's instructions.
5. Store the prismatic tube in such a way that it is protected from contamination until the next use.
6. Carefully re-attach the prismatic tube and allow it to lock in place (see section on "Refitting the prismatic tube [→ 17]").

3.2 CEREC Bluecam with sapphire glass

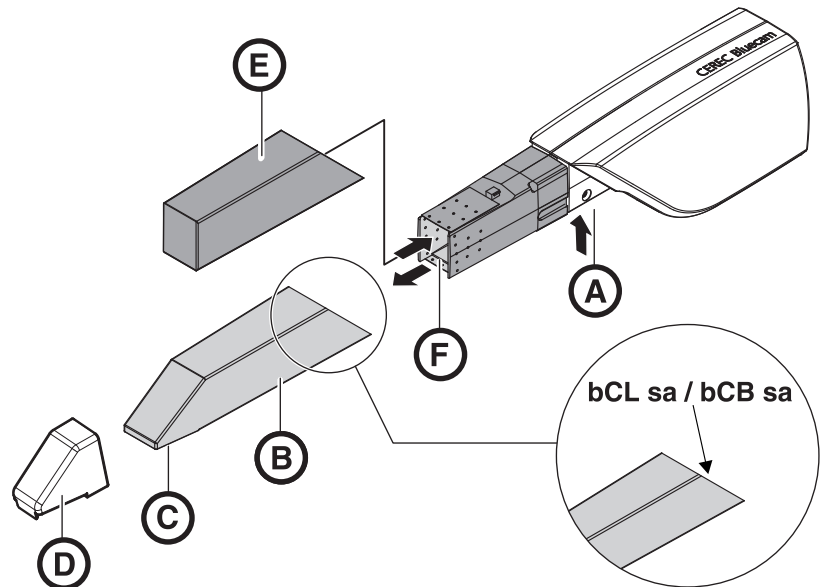
3.2.1 General information

The CEREC Bluecam is a very sensitive optical device and must therefore be handled with the **utmost care**. Protect the front lens or prism against scratches and clean them with a lint-free cloth and ethanol (commercially available cleaning alcohol) if any haze is noted during the acquisition.

The following methods are available for cleaning the prismatic tube:

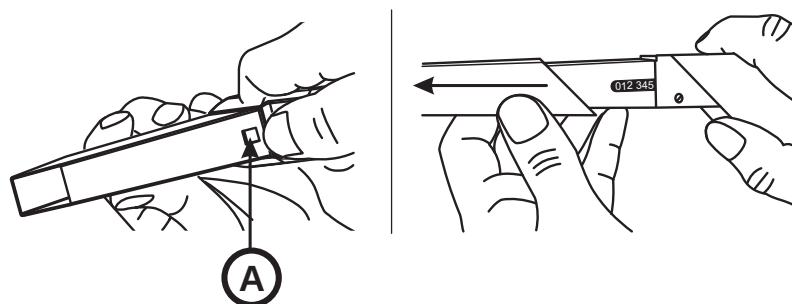
- Wipe disinfection for the camera and prismatic tube [→ 23]
- High-level disinfection of the prismatic tube without sapphire glass [→ 18]

3.2.2 Components of the CEREC Bluecam with sapphire glass

A	Press detent to release	D	Camera support 6 pcs, Order No. 5945360
B	Prismatic tube	E	Protective cap
C	Prism with sapphire glass (bCL sa / bCB sa)	F	Front lens

3.2.3 Removing the prismatic tube

A	Detent
---	--------

1. Press the prismatic tube against the camera body.
2. Press detent A.

NOTE

Risk of damaging the front lens or prism.

➤ Push the prismatic tube straight toward the front, do not tilt it.

3. Pull off the prismatic tube.

3.2.4 Refitting the prismatic tube

NOTE

Do not use CEREC 2 / CEREC 3 prismatic tubes.

NOTE

Risk of damaging the front lens or prism.

➤ The prismatic tube must not touch the front lens.
➤ Push the prismatic tube straight toward the camera body; do not tilt it.

- Carefully refit the prismatic tube until it locks in place.

3.2.5 Pre-cleaning of the prismatic tube

Cleaning process immediately after using the CEREC Bluecam: Wipe the camera with the appropriate prismatic tube so that the superficial impurities cannot harden and stick to the surface of the prismatic tube. Remove any visible contaminations using cotton wool gauze, moistened with tap water and a neutral soap, or lint-free cloth soaked with one of the disinfectant agents specified in the "Cleaning agents and disinfectants [→ 3]" section. Wipe down the prismatic tube afterwards again with the absorbent cotton gauze dipped in drinking water. Next dry the prismatic tube using a lint-free cloth.

3.2.6 Wipe disinfection for the camera and prismatic tube

NOTE

Do not spray the CEREC Bluecam with or immerse it in cleaning agents or disinfectants!

1. Remove the camera support (if present).
2. Carry out a preliminary cleaning process (see "Pre-cleaning of the prismatic tube [-> 22]").
3. Use a new absorbent cotton gauze which is moistened with one of the agents listed in the "Cleaning agents and disinfectants [-> 3]" section. First wipe the camera enclosure and then the prismatic tube. The camera enclosure and prismatic tube are disinfected through this step.
4. Wipe down the prismatic tube afterwards again with the absorbent cotton gauze dipped in drinking water.
5. Next dry the prismatic tube using a lint-free cloth.

3.2.7 High-level disinfection of the prismatic tube with sapphire glass

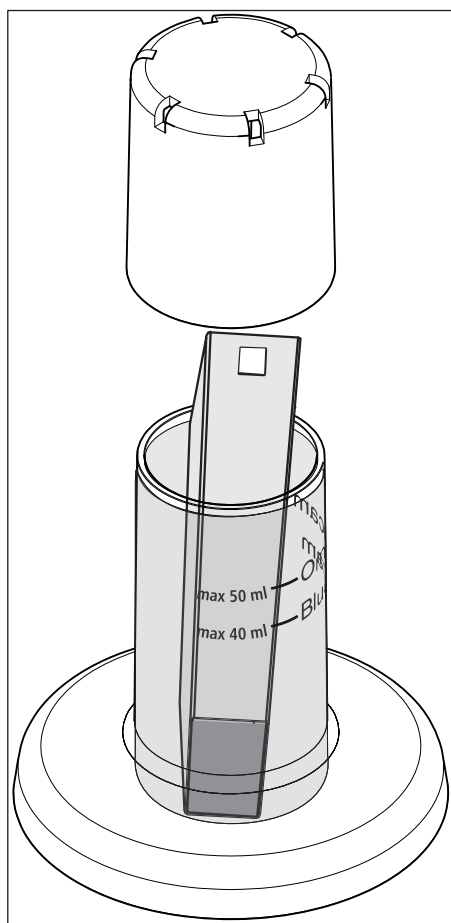
The complete process for high-level disinfection (HLD) is as follows – assuming the CIDEX[®] OPA disinfectant is available in your country via your dealer:

A HLD set to support the HLD processes can be ordered from Dentsply Sirona with REF 63 46 907.

1. Remove the camera support (if present).
2. Carry out a preliminary cleaning process (see "Pre-cleaning of the prismatic tube [→ 22]").
3. Always use personal protective equipment when implementing the high-level disinfection if not before.
4. Remove the prismatic tube from the CEREC Bluecam (see "Removing the prismatic tube [→ 22]").
5. Place the protective cap on the camera head and place the camera in the camera storage unit.
6. Use the following disinfectants for the high-level disinfection: CIDEX[®] OPA.

CAUTION

Observe the safety instructions from the disinfectant manufacturer.

7. Carefully fill the HLD container up to the 40 ml mark ("Bluecam" mark). A funnel can be used to help pour in the liquid. Do not spill the disinfectant. Should you spill any disinfectant, follow the safety instructions from the disinfectant manufacturer.

NOTE

Only place the prismatic tube into the container and not the entire camera.

Ensure that the prismatic tube is placed upright in the disinfectant (with the prism facing down) so that the solution cannot penetrate the inside of the prismatic tube.

8. Insert the prismatic tube in the correct position (with the prism facing down).
9. Place the cover on the container and leave the prismatic tube in the disinfectant for at least 12 minutes (CIDEX[®] OPA). The prismatic tube will not be damaged if it remains in the disinfectant for longer than the required time, but removing the prismatic tube shortly after the required time is recommended.
10. Remove the prismatic tube from the container.
11. Carefully rinse the prismatic tube with tap water for at least 30 seconds. Ensure that no water is able to penetrate into the interior of the prismatic tube.
12. Wipe down the prismatic tube. Use cotton wool gauze to do so, moistened with tap water.
13. Dry the prismatic tube using a soft, lint-free cloth.
14. Store the prismatic tube in such a way that it is protected from contamination until the next use.
15. Remove the protective cap from the camera head before use.

16. Carefully re-attach the prismatic tube and allow it to lock in place.
Refer to the section on "Refitting the prismatic tube [→ 22]".

We reserve the right to make any alterations which may be required due to technical improvements.

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© Sirona Dental Systems GmbH
Sirona Dental Systems GmbH
Sirona Dental, Inc.